

King Vision® has been designed with the intention to make a revolutionary series of high performance portable video laryngoscopes. The King Vision® combines the convenience of a durable reusable, video system with a disposable blade at a price that's affordable for all of your intubations, difficult and routine.

Learn more at www.airwayelearning.com

BETTER PATIENT CARE

The King Vision® is ergonomically designed to provide minimal lifting of soft tissue and impact on teeth.



PRODUCT SPECIFICATIONS

REUSABLE DISPLAY/ VIDEO ADAPTER				DISPOSABLE BLADES	
Batteries	3 AAA (Alkaline recommended)	Material	Polycarbonate / ABS anti-reflective coating on display window	Anti-fog	Coating on distal lens
Battery Life	> 90 min. (Battery status indicator flashes red when batteries need to be changed.)	Camera Chip	CMOS	AP Height of Blade	13 mm (non-channeled); 18 mm (channeled)
Computerized Power Management System	Auto shut off, Auto white balancing	Camera Resolution	640 x 480 VGA	Blade Length	17 cm
Video Aspect Ratio	4:3	Light Source	White LED		
Video Port	RCA connection to monitor with cable			Blade Width	26 mm (standard); 29 mm (channeled); 16 mm at distal tip
Video Refresh Rate	30 frames per second			Disposable Material	polycarbonate/TPE
Video Resolution	320 x 240 (QVGA)			ET Size	6.0 mm - 8.0 mm (Channeled Blade)
Video Screen	OLED Display				
Video Screen Size	6.1 cm / 2.4" diagonal				

ORDERING SPECIFICATIONS

ITEM NO.	DESCRIPTION	QUANTITY
KVLAB3	Standard aBlade, Size 3	20/CS
KVLAB3C	Channeled aBlade, Size 3	20/CS
KVLVA34	aBlade Video Adapter 3-4	1/CS
KVLABKIT3	1 Vision Display, 1 Video Adapter, Case, AAA Batteries, 4 Channeled and 2 Standard aBlade	1/CS
KVIS01VA34	1 Display and 1 Video Adapter 3-4	1/CS
KVL03	Standard Blade, Size 3	10/CS
KVL03C	Channeled Blade, Size 3	10/CS
KVIS01	Reusable Display	1/CS
KVLKIT3	Vision Display, Case, AAA Batteries, 3 Channeled and 1 Standard Size 3 Vision Blades	1/CS
KVCABL	Custom Video-out Cable (2.75 m / 9')	1/CS
KVIS01ADCASE	aBlade case	1/CS

CE US: Rx only

Sales Company
Ambu A/S
Baltorpbakken 13
DK - 2750 Ballerup
Denmark
Tel. +45 7225 2000
Fax +45 7225 2053
www.ambu.com

Manufacturer
King Systems
15011 Herriman Boulevard
Noblesville, IN 46060
Tel. 317 776 6823
Tel. 800 642 5464
Fax 317 776 6827
www.kingsystems.com



The growing challenge in healthcare around the world fuels our work to release the power of single-use devices and transform ideas into breakthrough solutions that work for life.

496 9901 01 • Voz • 2014/05 • Ambu A/S



King Vision® aBlade video laryngoscope

Leading the way in visualisation

Portable. Affordable. Durable.

King Vision® video laryngoscope


- Advanced imaging technology and unparalleled ergonomic design. Visionary performance in the palm of your hand.
- Unique flexibility from the choice of two blade types - standard and channelled - depending on the user technique and complexity of the airway. Choose what is right for you.
- King Vision® is designed to visualise the vocal cords indirectly which can help minimize the need for tissue manipulation and improve the ability to see anterior airways.

PORTABLE

- King Vision® is light weight, self-contained, battery operated and water resistant.
- King Vision® can be used continuously for more than 90 minutes making it highly portable.
- Blades are individually packaged. Each blade is disposable, eliminating concerns about cross contamination.
- King Vision® has a High Intensity OLED display that creates clear image viewing in 160° panoramic field which is also beneficial in teaching situations.

AFFORDABLE

- Disposable aBlades allow economical use of the King Vision® video laryngoscope for all intubations. With the addition of the aBlade configuration you have an excellent solution to the problem of limited availability of advanced imaging devices in the OR.
- With the King Vision® you can improve patient care by having high performance visualisation capabilities for every intubation, at an affordable price.

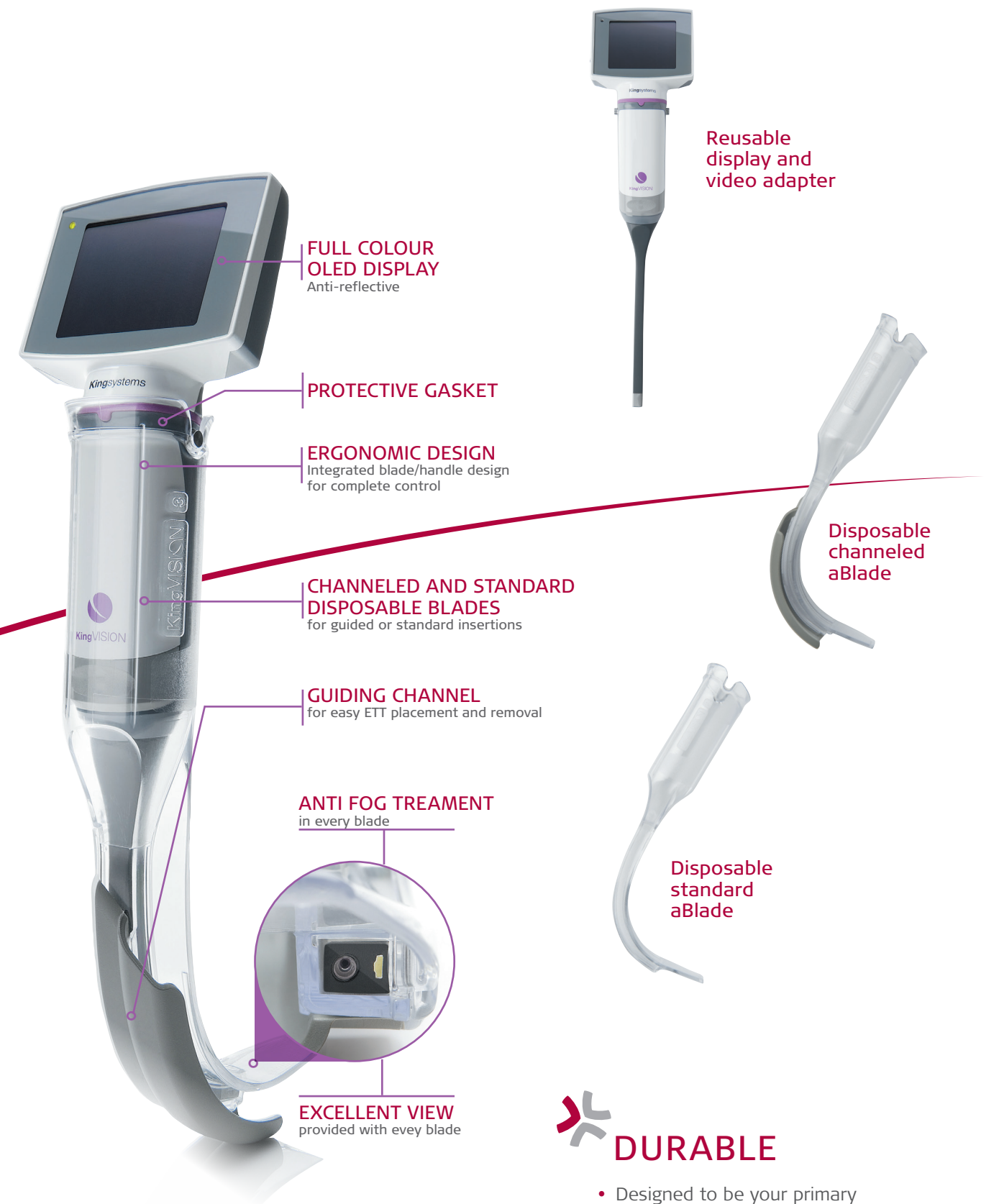

DURABLE

- Designed to be your primary tool for all intubations.
- Video output capability is compatible with external monitors and recording devices.
- The robust, full-colour, non-glare display can resist repeated cleaning and normal use wear and tear. The video adapter is made from durable materials that can withstand repeated use.