

VivaSight Max 7" Monitor



Key Benefits

- Enabling real-time high resolution video image
- Enabling continuous visualisation and monitoring throughout the Single Lung Ventilation procedure

VivaSight Max 7" Monitor

Continuous monitoring during Single Lung Ventilation

Intubation and airway management during isolation procedures present a series of challenges for even the most experienced team. Visualisation during Single Lung Ventilation (SLV) increases patient safety and is increasingly recommended as good clinical practice.

Single Lung Ventilation made simple and safe

VivaSight-DL and VivaSight-SL are the new generation of endotracheal tubes with integrated high-resolution camera. They make placement of double lumen tube and bronchial blocker fast and effective and the continuous visual monitoring

throughout the procedure makes single lung ventilation safer as incorrect position of the tube or bronchial blocker is easily detected and corrected.

The real-time high resolution image from the VivaSight-DL and VivaSight-SL is displayed on the medical grade LCD VivaSight Max 7" monitor. The monitor is 18.5 cm x 12.2 cm x 3.2 cm in size.

The full Ambu portfolio for SLV also includes VivaSight-DL (DLT) with integrated camera, VivaSight-SL (ETT) with integrated camera, VivaSight-EB bronchial blocker, and the single-use bronchoscope aScope™ 3 Slim.

Specifications

Size
7", 18.5 cm x 12.2 cm x 3.2 cm
Image
Resolution: 800 x 480 320 x 240 CIF high resolution image (76.800 pixels)
Video
Video Input Mode : NTSC Video Input Signal: CVBS (Composite Video Baseband Signal)
Monitor
Type: TFT LCD Response time: max. 30 ms 12Vdc 1.0A from an AC/DC adaptor
Power supply input
100-240 VAC 50-60 Hz, 0.3A
Extension cable
(Stereo Jack 4-pin) for power and signal transmission, 12V medical grade power supply


Ordering specifications

Description	Item No.	Quantity /box
VivaSight Max 7" LCD monitor	SA00014	1
Video Extension Cable for VivaSight Max monitor	SA00030	10
VivaSight Direct RCA cable kit for external monitors	SA00013	1
Power supply for VivaSight Max monitor (12V Medical)	EL0001700	1



VivaSight Max 7" Monitor
Compatible with VivaSight-DL
and VivaSight-SL

+

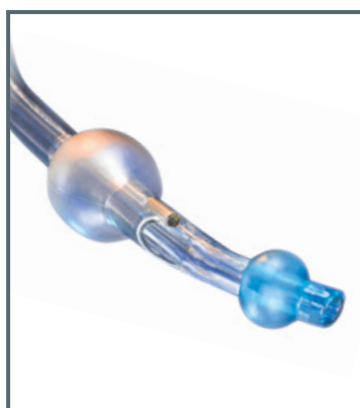

VivaSight-DL
Double Lumen Tube with integrated
camera in tracheal lumen

OR



VivaSight-SL
Single lumen endotracheal tube with
integrated camera at tip of tube.

Manufactured
by **ETView Ltd.**
Catom 2 Street
Misgav Business Park
M.P. Misgav 2017900
Israel

EU Representative: **Qsite EU**
Gerrit van der Veenstraat 84HS
1077 EL Amsterdam
The Netherlands
Tel: +31 (0) 20 811 0550
Fax: +31 (0) 84 221 3142
qsiteeu@qsite-med.com

VivaSight is distributed by **Ambu A/S**
Ambu A/S
Baltorpbakken 13
DK-2750 Ballerup
Denmark

Rx ONLY

CE 0483